

Ancor Pro™

Temporary Orthodontic
Anchorage System

The most innovative treatment option takes a new turn.



Yours.

Control and confidence at every turn.

Establishing and maintaining effective anchorage is key to gaining control over both the quality of the results and the duration of many of your orthodontic treatments.

The Ancor Pro Orthodontic Anchorage System by Ortho Organizers allows you to have greater control of patient outcomes by providing increased stabilization compared to traditional stabilization techniques.

With the turn of a screw, the Ancor Pro Anchor delivers absolute anchorage, and decreased reliance on patient compliance.

Superior versatility, simplicity and ease of use – by design.

Developed by a team with more than 25 years of implant industry experience, the Ancor Pro Orthodontic Anchorage System has been specially designed to provide anchorage in a wide range of clinical applications, while offering simplicity, ease of use and patient comfort.

Self-tapping & self-drilling mechanics – sharp screw tips and threads allow for easy placement and removal chairside without general anesthesia.*

Multi-functional single head – supragingival head features an upper button and undercut for the attachment of elastic chain, closed loop coil springs and other auxiliary orthodontic devices. Tapered conical top features a lumen to accept wires up to .022" in diameter. This innovative multi-functionality means greater procedural versatility and streamlined inventory control.


Variety of popular sizes – available in three lengths (6mm, 8mm, 10mm), and 1.6mm diameter; made from Grade 5 Titanium; providing maximum strength for ultimate performance.



*Local anesthesia is suitable, in most cases

The absolute anchorage system.

Complete system – the Ancor Pro Orthodontic Anchorage System includes an autoclavable, fully loaded kit that allows for easy storage of all components while not in use.

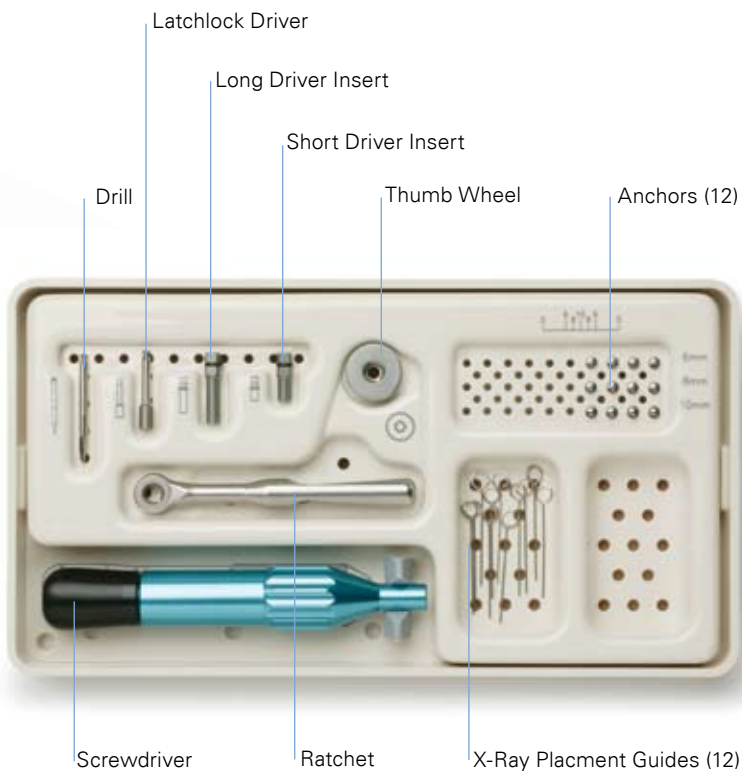

Ancor Pro Orthodontic Anchorage System provides a new means for increased stabilization in a variety of treatments. Moreover, it is backed by the rigorous quality control you have come to expect from Ortho Organizers.

Rely on Ancor Pro Orthodontic Anchorage System with confidence for the following treatments:

- Distalization of molars in the mandible and maxilla
- Mesialization of individual teeth or entire groups of teeth
- “En-masse” retraction
- Class II sliding mechanics
- Asymmetric Class II
- Space closure in a Class I occlusion
- Asymmetric occlusal plane

Due to the increased load tolerance and simplified biomechanics of the Ancor Pro Orthodontic Anchorage System, you can now establish immediate and absolute anchorage in locations previously impossible without invasive procedures such as orthognathic surgery and/or the extraction of teeth. This is also an alternative to the need for traditional devices such as headgear and other extraoral appliances.

Less treatment time. Less reliance on adherence. More support.

With its advanced design, innovative features and proven benefits, the Ancor Pro Orthodontic Anchorage System will increase patient comfort and help shorten treatment time.

Ancor Pro Orthodontic Anchorage System offers you the added benefit of unparalleled support by Ortho Organizers – responsive customer service; knowledgeable and experienced staff; and comprehensive educational offerings.

Catalog Number	Part Name
500-000	Ortho Anchor Kit Anchors (12) Sterilization Kit Screwdriver Long Driver Insert Short Driver Insert Latchlock Driver Thumb Wheel Drill X-Ray Placement Guides (12) Comfort Caps (12)**
500-006	6 mm anchor
500-008	8 mm anchor
500-010	10 mm anchor
500-020	Drill
500-021	Latchlock Driver
500-022	Short Driver Insert
500-023	Long Driver Insert
500-024	Thumb Wheel
500-025	Ratchet
500-026	X-Ray Placement Guides
500-027	Screwdriver
500-028	Comfort Caps**

For more information on our products and educational offerings, please contact us:
 In the U.S.: 800.547.2000
 Outside the U.S.: +1.760.448.8600
 Canada: CERUM. 800.661.9567
 To fax an order: 1.800.888.7244
 To email an order: usasales@OrthoOrganizers.com

Ancor Pro™
Temporary Orthodontic
 Anchorage System



1822 Aston Avenue
 Carlsbad, CA 92008
 USA
www.OrthoOrganizers.com