
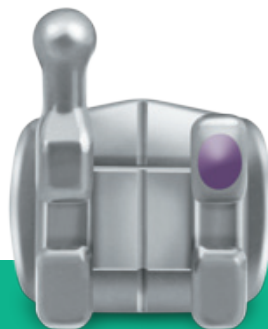


Maestro™

Bracket System

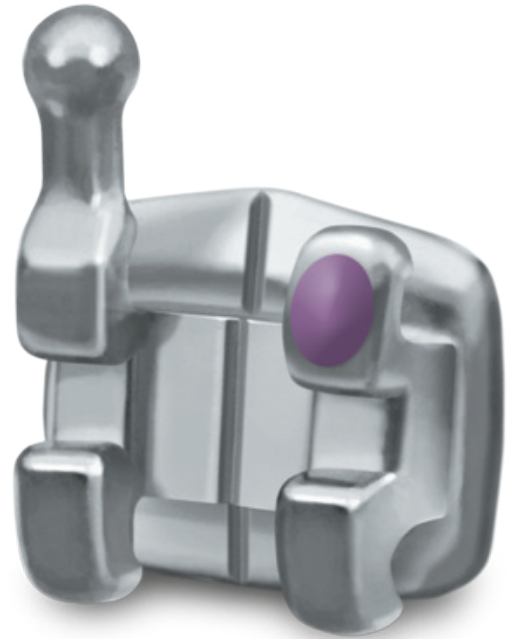


Discover the Precision, Comfort, and Control of the MIM-Manufactured Maestro Bracket

The *Maestro*™ Low-Profile Bracket System harmonizes the enduring precepts of classic orthodontic design with state-of-the-art engineering and production practices. The result is a complete orthodontic appliance system with corresponding buccal tubes—contemporary in form and function—that achieves the optimum in metal aesthetics and comfort that patients seek without sacrificing the precision and control that practitioners need. A finely tuned appliance design for predictable, efficient treatment and discerning finishes.


Metal Injection Molding (MIM) is advanced technology for manufacturing complex metal products. Fine metal powders and wax binders are injected into intricate molds. After the wax binders are removed, the brackets are sintered to form a high-density orthodontic

appliance with extremely close tolerances that accurately replicates the intricacies of the original mold. MIM translates into rounded contours, and complex structures that machined brackets cannot duplicate providing corrosion resistance, strength, and proven performance.

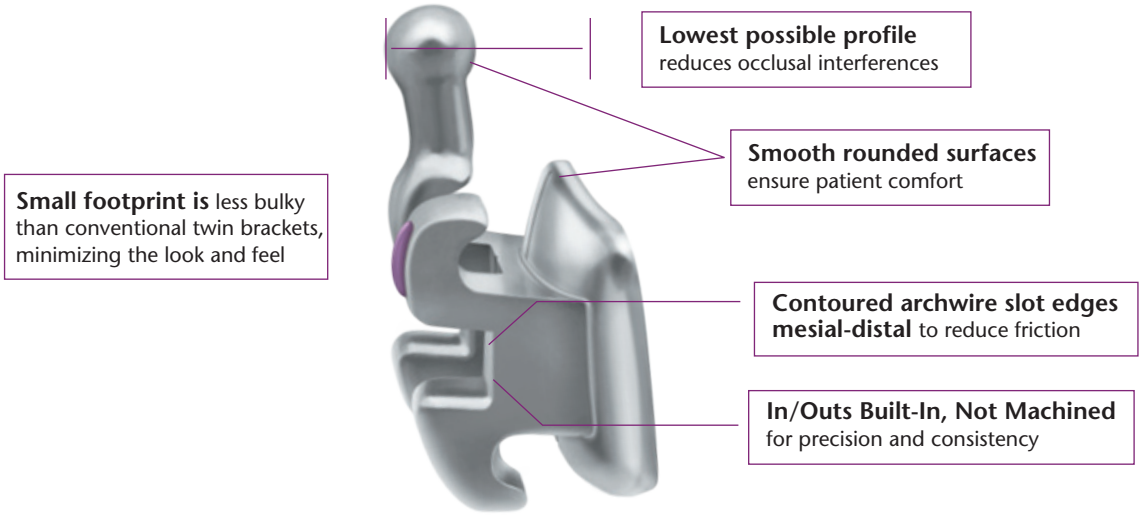

THE MAESTRO BRACKET OFFERS MULTIPLE VISUAL CUES TO ENHANCE PRECISE BRACKET PLACEMENT

The geometries of the *Maestro* Bracket are specifically designed to provide signposts which assist you in positioning brackets in their ideal location.

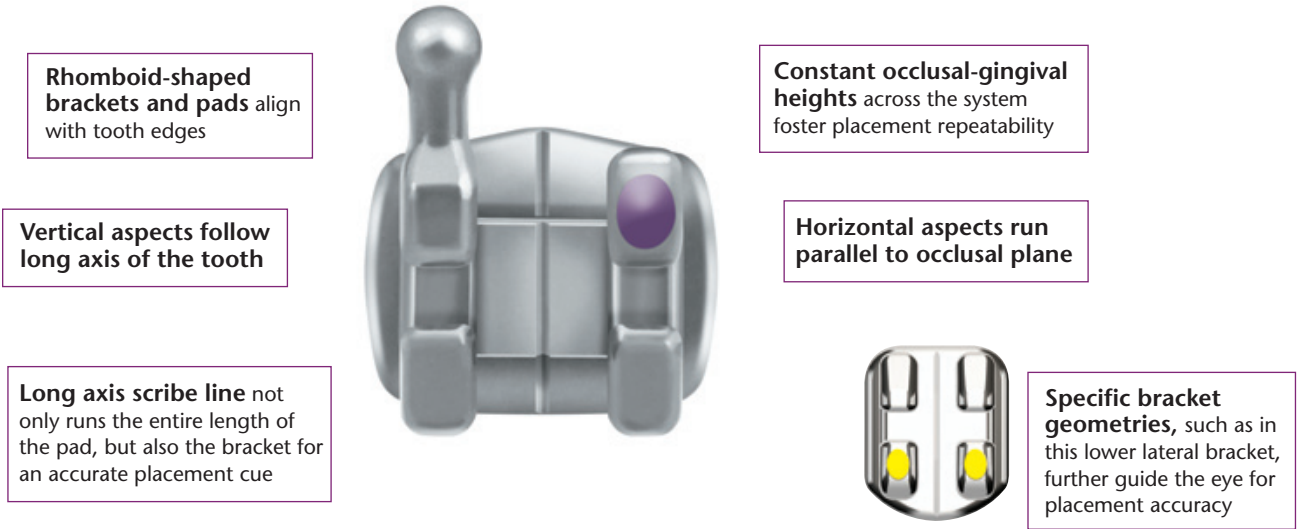

Maestro: A Complete System

Maestro Buccal Tubes complete the appliance system. The 80 gauge mesh pad, compound contoured base, and defined buccal indent, provides the practitioner with the precise control needed throughout all phases of treatment. The large funneled entrance, side-grip areas, and color-coding make for easier identification, placement, positioning, plus wire insertion. The patient will have one of the most comfortable experiences with its smooth shape, low profile, rounded edges, and lowered hook. The *Maestro* Buccal Tube orchestrates all components of the system seamlessly, providing reduced chair time and increased satisfaction for both practitioner and patient.



Maestro Buccal Tubes—bondable and weldable—are engineered to provide a complete and seamless appliance system

Funneled slot opening makes wire engagement easier without adding bulk

IN-DEPTH RESEARCH ON TOOTH MORPHOLOGY INFORMS SUPERIOR ENGINEERING FINELY TUNED DESIGN FOR DETAILING FINELY TUNED FINISHES



Anatomically contoured base and pad nestle each tooth for positioning that *feels right*



80 gauge micro-etched bondable mesh pad offers proven bond adhesion



Small profile is balanced with **deep tie-wing undercuts** that ease secure ligation



Torque-in-base provides level slot lineup, and helps avoid occlusal interferences



Optional **gingivally offset brackets, on large bonding bases** of bicuspid, minimize bond failures



Optimized mesial-distal width offers superb rotational control

The *Maestro* Bracket is available in the Roth prescription, and the McLaughlin, Bennett, and Trevisi prescriptions, in .018" and .022" lumen sizes. Ball hooks are optional on cuspids and bicuspid brackets; and each bracket is color-coded with an ID mark for easy identification.

MAESTRO MINI-TWIN® BRACKETS












ALL SETS (5x5)*

Rx	DESCRIPTION	QTY OF SETS	TOTAL BRACKETS	ITEM NUMBERS	
				.018	.022
Roth	Set with no hooks	1	20	711-297	711-397
	Set with no hooks, quick pac	1	20	711-297Q	711-397Q
	Set with no hooks	5	100	711-297-5	711-397-5
	Set with no hooks	10	200	711-297-10	711-397-10
	Set with no hooks & offset pads on lower 4's & 5's	1	20	711-297O	711-397O
	Set with no hooks & offset pads on lower 4's & 5's, quick pac	1	20	711-297OQ	711-397OQ
	Set with hooks on 3's	1	20	711-298	711-398
	Set with hooks on 3's, quick pac	1	20	711-298Q	711-398Q
	Set with hooks on 3's	5	100	711-298-5	711-398-5
	Set with hooks on 3's	10	200	711-298-10	711-398-10
	Set with hooks on 3's, 4's, 5's	1	20	711-299	711-399
	Set with hooks on 3's, 4's, 5's, quick pac	1	20	711-299Q	711-399Q
	Set with hooks on 3's, 4's, 5's	5	100	711-299-5	711-399-5
	Set with hooks on 3's, 4's, 5's	10	200	711-299-10	711-399-10
	Set with hooks on 3's, 4's, 5's & offset pads on lower 4's & 5's	1	20	711-299O	711-399O
	Set with hooks on 3's, 4's, 5's & offset pads on lower 4's & 5's, quick pac	1	20	711-299OQ	711-399OQ
	Set with hooks on 3's, 4's, 5's & offset pads	5	100	711-299O-5	711-399O-5
	Set with hooks on 3's, 4's, 5's & offset pads	10	200	711-299O-10	711-399O-10

Rx	DESCRIPTION	QTY OF SETS	TOTAL BRACKETS	ITEM NUMBERS	
				.018	.022
McLaughlin, Bennett, and Trevisi	Set with no hooks	1	20	711-297MC	711-397MC
	Set with no hooks, quick pac	1	20	711-297MCQ	711-397MCQ
	Set with no hooks	5	100	711-297MC-5	711-397MC-5
	Set with no hooks	10	200	711-297MC-10	711-397MC-10
	Set with no hooks & offset pads on lower 4's & 5's	1	20	711-297MCO	711-397MCO
	Set with no hooks & offset pads on lower 4's & 5's, quick pac	1	20	711-297MCOQ	711-397MCOQ
	Set with hooks on 3's	1	20	711-298MC	711-398MC
	Set with hooks on 3's, quick pac	1	20	711-298MCQ	711-398MCQ
	Set with hooks on 3's	5	100	711-298MC-5	711-398MC-5
	Set with hooks on 3's	10	200	711-298MC-10	711-398MC-10
	Set with hooks on 3's, 4's, 5's	1	20	711-299MC	711-399MC
	Set with hooks on 3's, 4's, 5's, quick pac	1	20	711-299MCQ	711-399MCQ
	Set with hooks on 3's, 4's, 5's	5	100	711-299MC-5	711-399MC-5
	Set with hooks on 3's, 4's, 5's	10	200	711-299MC-10	711-399MC-10
	Set with hooks on 3's, 4's, 5's & offset pads on lower 4's & 5's	1	20	711-299MCO	711-399MCO
	Set with hooks on 3's, 4's, 5's & offset pads on lower 4's & 5's, quick pac	1	20	711-299MCOQ	711-399MCOQ

* Ortho Organizers® prescription does not imply endorsement by the doctor(s).

MAESTRO MINI-TWIN BRACKETS
ALL PRESCRIPTIONS*

MAXILLARY		DESC.	COLOR-CODE	TOOTH #	TORQUE	ANG.	ANTI-ROTATION	ITEM NUMBERS			
								.018 R	.018 L	.022 R	.022 L
Centrals		—	Black	U1	+12°	+5°	0°	711-203	711-204	711-303	711-304
		—	Black	U1	+17°	+4°	0°	711-205	711-206	711-305	711-306
Laterals		—	Pink	U2	+8°	+9°	0°	711-221	711-222	711-321	711-322
		—	Pink	U2	+10°	+8°	0°	711-223	711-224	711-323	711-324
Cuspids		—	Green	U3	-2°	+13°	+4°M	711-251	711-252	711-351	711-352
		—	Green	U3	-7°	+8°	0°	711-255	711-256	711-355	711-356
		Hook	Green	U3	-2°	+13°	+4°M	711-251HK	711-252HK	711-351HK	711-352HK
		Hook	Green	U3	-7°	+8°	0°	711-255HK	711-256HK	711-355HK	711-356HK
Bicuspids		—	Purple	U4&5's	-7°	0°	+2°D	711-259	711-260	711-359	711-360
		—	Purple	U4&5's	-7°	0°	0°	711-257	711-258	711-357	711-358
		Hook	Purple	U4&5's	-7°	0°	+2°D	711-259HK	711-260HK	711-359HK	711-360HK
		Hook	Purple	U4&5's	-7°	0°	0°	711-257HK	711-258HK	711-357HK	711-358HK
		Offset pad	Purple	U4&5's	-7°	0°	+2°D	711-479	711-480	711-579	711-580
		Offset pad	Purple	U4&5's	-7°	0°	0°	711-457	711-458	711-557	711-558
		Offset pad w/hook	Purple	U4&5's	-7°	0°	+2°D	711-479HK	711-480HK	711-579HK	711-580HK
		Offset pad w/hook	Purple	U4&5's	-7°	0°	0°	711-457HK	711-458HK	711-557HK	711-558HK
MANDIBULAR		DESC.	COLOR-CODE	TOOTH #	TORQUE	ANG.	ANTI-ROTATION	ITEM NUMBERS			
								.018 R	.018 L	.022 R	.022 L
Anteriors		—	Yellow	L1&2's	-1°	0°	0°	711-233	711-233	711-333	711-333
		—	Yellow	L1&2's	-6°	0°	0°	711-235	711-235	711-335	711-335
Cuspids		—	Blue	L3	-11°	+7°	+2°M	711-269	711-270	711-369	711-370
		—	Blue	L3	-6°	+3°	0°	711-271	711-272	711-371	711-372
		Hook	Blue	L3	-11°	+7°	+2°M	711-269HK	711-270HK	711-369HK	711-370HK
		Hook	Blue	L3	-6°	+3°	0°	711-271HK	711-272HK	711-371HK	711-372HK
1st Bicuspids		—	White	L4	-17°	0°	+4°D	711-289	711-290	711-389	711-390
		—	Gray	L4	-12°	+2°	0°	711-287	711-288	711-387	711-388
		Hook	White	L4	-17°	0°	+4°D	711-289HK	711-290HK	711-389HK	711-390HK
		Hook	Gray	L4	-12°	+2°	0°	711-287HK	711-288HK	711-387HK	711-388HK
		Offset pad	White	L4	-17°	0°	+4°D	711-279	711-280	711-379	711-380
		Offset pad	Gray	L4	-12°	+2°	0°	711-477	711-478	711-577	711-578
		Offset pad w/hook	White	L4	-17°	0°	+4°D	711-279HK	711-280HK	711-379HK	711-380HK
		Offset pad w/hook	Gray	L4	-12°	+2°	0°	711-477HK	711-478HK	711-577HK	711-578HK
2nd Bicuspids		—	Red	L5	-22°	0°	+4°D	711-293	711-294	711-393	711-394
		—	White	L5	-17°	+2°	0°	711-291	711-292	711-391	711-392
		Hook	Red	L5	-22°	0°	+4°D	711-293HK	711-294HK	711-393HK	711-394HK
		Hook	White	L5	-17°	+2°	0°	711-291HK	711-292HK	711-391HK	711-392HK
		Offset pad	Red	L5	-22°	0°	+4°D	711-473	711-474	711-573	711-574
		Offset pad	White	L5	-17°	+2°	0°	711-471	711-472	711-571	711-572
		Offset pad w/hook	Red	L5	-22°	0°	+4°D	711-473HK	711-474HK	711-573HK	711-574HK
		Offset pad w/hook	White	L5	-17°	+2°	0°	711-471HK	711-472HK	711-571HK	711-572HK

* Ortho Organizers prescription does not imply endorsement by the doctor(s).
 Bracket images are not representative of actual size, or ratio to other brackets.

MAESTRO MINI-TWIN BRACKETS

ALL PRESCRIPTIONS*

MAXILLARY	CENTRALS			LATERALS			CUSPIDS			1ST BICUSPIDS		2ND BICUSPIDS			
	TRO.	ANG.	ANTI-ROTATION	TRO.	ANG.	ANTI-ROTATION	TRO.	ANG.	ANTI-ROTATION	TRO.	ANG.	ANTI-ROTATION	TRO.	ANG.	ANTI-ROTATION
Roth	+12°	+5°	0°	+8°	+9°	0°	-2°	+13°	+4°M	-7°	0°	+2°D	-7°	0°	+2°D
McLaughlin, Bennett, and Trevisi	+17°	+4°	0°	+10°	+8°	0°	-7°	+8°	0°	-7°	0°	0°	-7°	0°	0°

MANDIBULAR	CENTRALS			LATERALS			CUSPIDS			1ST BICUSPIDS		2ND BICUSPIDS			
	TRO.	ANG.	ANTI-ROTATION	TRO.	ANG.	ANTI-ROTATION	TRO.	ANG.	ANTI-ROTATION	TRO.	ANG.	ANTI-ROTATION	TRO.	ANG.	ANTI-ROTATION
Roth	-1°	0°	0°	-1°	0°	0°	-11°	+7°	+2°M	-17°	0°	+4°D	-22°	0°	+4°D
McLaughlin, Bennett, and Trevisi	-6°	0°	0°	-6°	0°	0°	-6°	+3°	0°	-12°	2°	0°	-17°	2°	0°

MAESTRO MINI-TWIN BRACKETS

ROTH PRESCRIPTIONS*

MAXILLARY	DESC.	COLOR-CODE	TOOTH #	TORQUE	ANG.	ANTI-ROTATION	ITEM NUMBERS				
							.018 R	.018 L	.022 R	.022 L	
Centrals		—	Black	U1	+12°	+5°	0°	711-203	711-204	711-303	711-304
Laterals		—	Pink	U2	+8°	+9°	0°	711-221	711-222	711-321	711-322
Cuspids		—	Green	U3	-2°	+13°	+4°M	711-251	711-252	711-351	711-352
		Hook	Green	U3	-2°	+13°	+4°M	711-251HK	711-252HK	711-351HK	711-352HK
Bicuspids		—	Purple	U4&5's	-7°	0°	+2°D	711-259	711-260	711-359	711-360
		Hook	Purple	U4&5's	-7°	0°	+2°D	711-259HK	711-260HK	711-359HK	711-360HK
		Offset pad	Purple	U4&5's	-7°	0°	+2°D	711-479	711-480	711-579	711-580
		Offset pad w/hook	Purple	U4&5's	-7°	0°	+2°D	711-479HK	711-480HK	711-579HK	711-580HK

MANDIBULAR	DESC.	COLOR-CODE	TOOTH #	TORQUE	ANG.	ANTI-ROTATION	ITEM NUMBERS				
							.018 R	.018 L	.022 R	.022 L	
Anteriors		—	Yellow	L1&2's	-1°	0°	0°	711-233	711-233	711-333	711-333
Cuspids		—	Blue	L3	-11°	+7°	+2°M	711-269	711-270	711-369	711-370
		Hook	Blue	L3	-11°	+7°	+2°M	711-269HK	711-270HK	711-369HK	711-370HK
1st Bicuspids		—	White	L4	-17°	0°	+4°D	711-289	711-290	711-389	711-390
		Hook	White	L4	-17°	0°	+4°D	711-289HK	711-290HK	711-389HK	711-390HK
		Offset pad	White	L4	-17°	0°	+4°D	711-279	711-280	711-379	711-380
		Offset pad w/hook	White	L4	-17°	0°	+4°D	711-279HK	711-280HK	711-379HK	711-380HK
2nd Bicuspids		—	Red	L5	-22°	0°	+4°D	711-293	711-294	711-393	711-394
		Hook	Red	L5	-22°	0°	+4°D	711-293HK	711-294HK	711-393HK	711-394HK
		Offset pad	Red	L5	-22°	0°	+4°D	711-473	711-474	711-573	711-574
		Offset pad w/hook	Red	L5	-22°	0°	+4°D	711-473HK	711-474HK	711-573HK	711-574HK

* Ortho Organizers prescription does not imply endorsement by the doctor(s).
Bracket images are not representative of actual size, or ratio to other brackets.

MCLAUGHLIN, BENNETT, AND TREVISI PRESCRIPTIONS*

MAXILLARY		DESC.	COLOR-CODE	TOOTH #	TORQUE	ANG.	ANTI-ROTATION	ITEM NUMBERS			
								.018 R	.018 L	.022 R	.022 L
Centrals		—	Black	U1	+17°	+4°	0°	711-205	711-206	711-305	711-306
Laterals		—	Pink	U2	+10°	+8°	0°	711-223	711-224	711-323	711-324
Cuspids		—	Green	U3	-7°	+8°	0°	711-255	711-256	711-355	711-356
		Hook	Green	U3	-7°	+8°	0°	711-255HK	711-256HK	711-355HK	711-356HK
Bicuspids		—	Purple	U4&5's	-7°	0°	0°	711-257	711-258	711-357	711-358
		Hook	Purple	U4&5's	-7°	0°	0°	711-257HK	711-258HK	711-357HK	711-358HK
		Offset pad	Purple	U4&5's	-7°	0°	0°	711-457	711-458	711-557	711-558
		Offset pad w/hook	Purple	U4&5's	-7°	0°	0°	711-457HK	711-458HK	711-557HK	711-558HK

MANDIBULAR								ITEM NUMBERS			
								.018 R	.018 L	.022 R	.022 L
Anteriors		—	Yellow	L1&2's	-6°	0°	0°	711-235	711-235	711-335	711-335
Cuspids		—	Blue	L3	-6°	+3°	0°	711-271	711-272	711-371	711-372
		Hook	Blue	L3	-6°	+3°	0°	711-271HK	711-272HK	711-371HK	711-372HK
1st Bicuspids		—	Gray	L4	-12°	+2°	0°	711-287	711-288	711-387	711-388
		Hook	Gray	L4	-12°	+2°	0°	711-287HK	711-288HK	711-387HK	711-388HK
		Offset pad	Gray	L4	-12°	+2°	0°	711-477	711-478	711-577	711-578
		Offset pad w/hook	Gray	L4	-12°	+2°	0°	711-477HK	711-478HK	711-577HK	711-578HK
2nd Bicuspids		—	White	L5	-17°	+2°	0°	711-291	711-292	711-391	711-392
		Hook	White	L5	-17°	+2°	0°	711-291HK	711-292HK	711-391HK	711-392HK
		Offset pad	White	L5	-17°	+2°	0°	711-471	711-472	711-571	711-572
		Offset pad w/hook	White	L5	-17°	+2°	0°	711-471HK	711-472HK	711-571HK	711-572HK

BRACKET IDENTIFICATION CHART

COLOR-CODED ID FOR EASY REFERENCE

● U5R ● U4R ● U3R ● U2R ● U1R ● U1L ● U2L ● U3L ● U4L ● U5L

● L4 & 5: -12°: GRAY ○ -17°: WHITE ● -22°: RED

● L4 & 5: -12°: GRAY ○ -17°: WHITE ● -22°: RED

● L5R' ● L4R' ● L3R ● L2R ● L1R ● L1L ● L2L ● L3L ● L4L' ● L5L'

MAXILLARY

U1: BLACK - DISTAL-GINGIVAL TIE WING
 U2: PINK - DISTAL-GINGIVAL TIE WING
 U3: GREEN - DISTAL-GINGIVAL TIE WING
 U4 & 5: PURPLE - GINGIVAL TIE WINGS = UNIVERSAL

MANDIBULAR

L1 & 2: YELLOW - GINGIVAL TIE WINGS = UNIVERSAL
 L3: BLUE - DISTAL-GINGIVAL TIE WING
 L4 & 5: SEE ABOVE

ROTH

LOWER 4'S: WHITE
 LOWER 5'S: RED





MCLAUGHLIN, BENNETT AND TREVISI

LOWER 4'S: GRAY
 LOWER 5'S: WHITE

* Ortho Organizers prescription does not imply endorsement by the doctor(s).
 Bracket images are not representative of actual size, or ratio to other brackets.





MAESTRO BUCCAL TUBES

MCLAUGHLIN, BENNETT, AND TREVISI PRESCRIPTIONS*

STD	DESC.	COLOR-CODE	TORQUE	ANG.	OFFSET	ITEM NUMBERS			
						WELDABLE	.018 DIRECT BOND	WELDABLE	.022 DIRECT BOND
	Upper Right	Purple	-14°	0°	+10°	902-203	902-203DB	902-303	902-303DB
	Upper Left	Purple	-14°	0°	+10°	902-204	902-204DB	902-304	902-304DB
	Lower Right	White	-20°	0°	0°	902-253	902-253DB	902-353	902-353DB
	Lower Left	White	-20°	0°	0°	902-254	902-254DB	902-354	902-354DB

MAESTRO BUCCAL TUBES

ROTH PRESCRIPTIONS*

STD	DESC.	COLOR-CODE	TORQUE	ANG.	OFFSET	ITEM NUMBERS			
						WELDABLE	.018 DIRECT BOND	WELDABLE	.022 DIRECT BOND
	Upper Right	Purple	-14°	0°	+14°	902-201	902-201DB	902-301	902-301DB
	Upper Left	Purple	-14°	0°	+14°	902-202	902-202DB	902-302	902-302DB
	Lower Right	Red	-25°	0°	+4°	902-251	902-251DB	902-351	902-351DB
	Lower Left	Red	-25°	0°	+4°	902-252	902-252DB	902-352	902-352DB

The *Maestro* Bracket System is manufactured in the United States, emphasizing *Ortho Organizers* dedication to providing world-class products to meet your every orthodontic need. As a full service partner, we offer solutions that diversify your treatment options and help expand your practice.

For more information on our products, orthodontic laboratory services, and educational offerings, please contact us:
 Customer Service: 888.851.0533
 In the U.S.: 800.547.2000
 Outside the U.S.: +(1) 760 448 8600
 Canada: CERUM 800.661.9567
 To fax an order: 800.888.7244
 To fax an order outside the U.S.: +(1) 760 448 8607
 To email an order: USASales@OrthoOrganizers.com
 To email an order outside the U.S.: IntlOrders@OrthoOrganizers.com

* *Ortho Organizers* prescription does not imply endorsement by the doctor(s).



1822 Aston Avenue
 Carlsbad, CA 92008
 USA
 OrthoOrganizers.com