



Dr. William J. Clark

TransForce²

Arch Developer Appliances

The next generation in expansion

TransForce² Arch Developers are available in two styles, Transverse and Sagittal. Each of these appliances can be used to treat all classes of malocclusion and are excellent for mixed dentition. Both styles mount lingually, making them very inconspicuous and require zero patient compliance. The *TransForce²* Arch Developer Appliances were developed in consultation with Dr. William J. Clark, the inventor of the Clark Twin Block[®].

The enhanced *TransForce²* Arch Developer Appliances offer the following advantages:

Ease of Insertion

Placement blades are designed to insert easily into molar bands with standard .036 x .072 lingual sheaths

Aesthetic Appearance

Lingual placement for a virtually invisible appearance

Size Identification

Color-coded for easy identification

Optimum Patient Comfort

Nitanium[®] Coil Springs enclosed in expansion modules provide slow continual force

No Laboratory Work Required


Select and order appliance by desired finish size

Less Inventory

Universal sizes for upper or lower arches – no need to carry separate inventory

Zero Patient Compliance

Appliance allows for up to 2mm expansion every 6-8 weeks with minimum maintenance

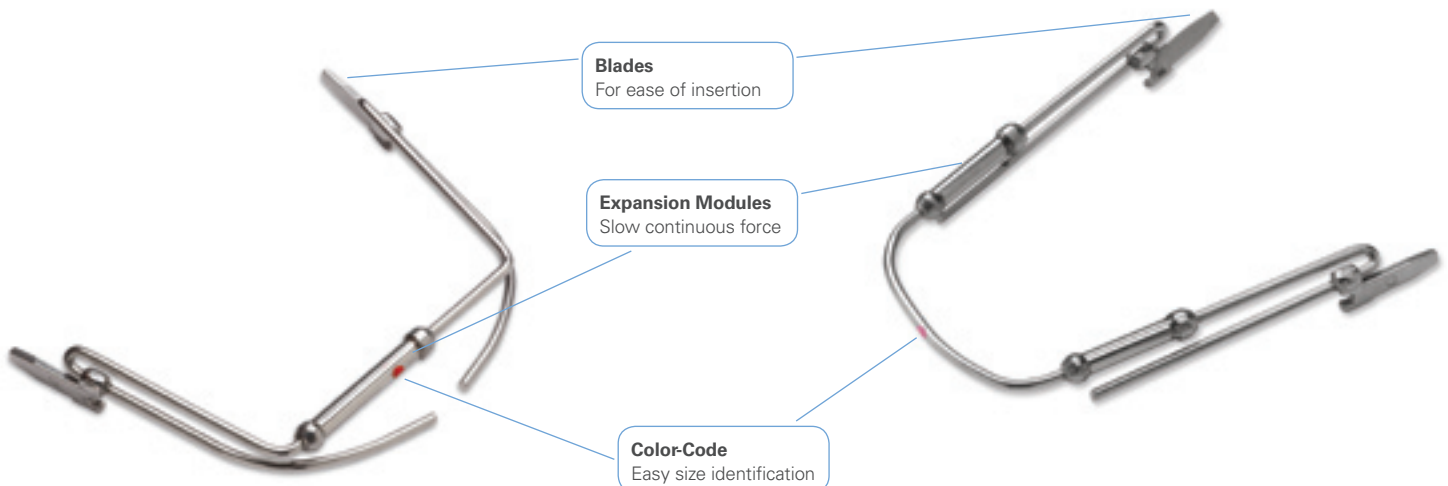

TransForce² Transverse Arch Developer

The Transverse Appliance has an expansion module to increase the inter-canine width (expands at the cuspids) and may be used in the upper or lower arch when expansion is required to accommodate crowding in the labial segments, or to correct arch width in constricted arches.

TransForce² Sagittal Arch Developer

The Sagittal Appliance is specifically designed for anterior arch development (creates space by advancing the 3x3 labially) and is often indicated for simultaneous use in both arches. The appliance operates on the slide principle and may be used unilaterally or bilaterally to extend arch length.





Virtually Invisible Assistance for Arch Development

Designed for ease of use, the *TransForce²* Arch Developer Appliances are preprogrammed to allow appliances to be placed without laboratory work. Both the Transverse and Sagittal styles may be placed using first molar bands with standard .036 x .072 lingual sheaths.

Each *TransForce²* Arch Developer Appliance may be easily adjusted per patient. In addition, each expander module contains *Nitanium[®]* Coil Springs generating approximately 200 grams of gentle biocompatible force.

The Transverse and Sagittal Appliances are available in multiple sizes, individually or a starter assortment kit. Best of all, the developers are universal, which allows for less inventory as they can be used for either upper or lower arch development.

TransForce² Transverse Arch Developer

Transverse Size	Catalog Number	Color-Code
26mm	424-526	 Red
28mm	424-528	 Green
30mm	424-530	 Purple
32mm	424-532	 Pink
Transverse Kit (Includes 1 of each size)	424-500	N/A

TransForce² Sagittal Arch Developer

Sagittal Size	Catalog Number	Color-Code
28mm	424-728	 Green
30mm	424-730	 Purple
32mm	424-732	 Pink
34mm	424-734	 Blue
36mm	424-736	 Black
38mm	424-738	 Yellow
40mm	424-740	 White
Sagittal Kit (Includes 1 of each size)	424-700	N/A

A *TransForce²* Arch Developer Planner may also be viewed on our Web site at www.OrthoOrganizers.com to assist with size selection.

For more information on our products, orthodontic laboratory services and educational offerings, please contact us:

Customer Service: 888.851.0533

In the U.S.: 800.547.2000

Outside the U.S.: +1.760.448.8600

Canada: CERUM. 800.661.9567

To fax an order: 800.888.7244

To fax an order outside the U.S.: +1.760.471.9549

To email an order: UsaSales@OrthoOrganizers.com



1822 Aston Avenue
Carlsbad, CA 92008
USA

www.OrthoOrganizers.com